

Direct-Q[®] 3, 5, 8 Water Purification Systems

Tap to pure and ultrapure water – with easy
and convenient dispense!



Tap to pure and ultrapure water – with easy and convenient dispense!

Your water purification needs

High quality water produced directly from tap water

Ultrapure water easily accessible wherever you need it in your lab

Compact all-in-one system for the most efficient use of your lab space

A choice of adequate storage volumes

Flow rates adapted to your daily usage

Instant flow rates to match your needs

High quality water to meet the requirements of your most critical applications

Easily accessible information on system operation

Simple, low-level self-maintenance

Our solution: the Direct-Q® range of water purification systems

Direct-Q® systems deliver **both pure and ultrapure water** directly from tap for use with a wide variety of applications in your lab.

With the **Direct-Q® range** of water purification systems, you benefit from a choice of ultrapure water dispensing possibilities. The innovative, space-saving **Remote dispenser** offers you water delivery solutions to best fit the way you work, with **easy and convenient remote delivery** up to two meters away from your water production unit.

The **Direct-Q®** range includes **versatile solutions** that can be bench-integrated or bench- / wall-installed.

Direct-Q® systems come either with an **integrated 6-liter reservoir** or an **external reservoir** (30 l or 60 l).

Systems in the Direct-Q® range provide **3, 5 or 8 liters of pure water per hour**, covering lab needs of 50,100 and 150 l of pure water per day.

Direct-Q® systems dispense up to 30 l of **ultrapure water per hour**.


Options such as a **UV lamp** and a range of **Application Pak point-of-use polishers** are available to fine-tune your ultrapure water.

The **user-friendly display** provides system status at a glance; the concise **Quick Reference Guide** is a handy guide for daily operation.

All-in-one SmartPak® cartridges enable easy and rapid replacement.

Direct-Q® Systems Water Purification Pathway

1. Booster Pump
2. Inlet Solenoid Valve
3. SmartPak® Part 1 (Pretreatment and RO Cartridge)
4. RO Reject Solenoid Valve
5. RO Reject Capillary
6. Check valve
7. RO Permeate Conductivity Cell
8. Reservoir (Built-in or separate, depending on model)
9. Distribution Pump
10. UV Lamp 185/254 nm (UV System)
11. SmartPak® Part 2 (Synthetic Activated Carbon & Ion Exchange Polisher Cartridge)
12. Product Resistivity Cell

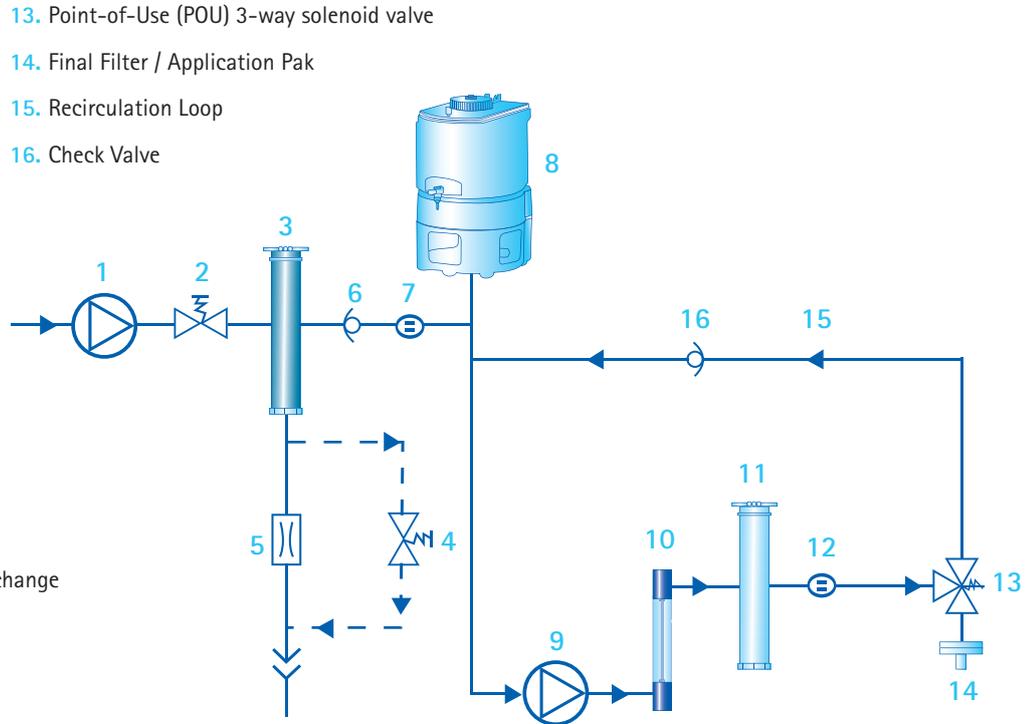

Choose the solution that's right for you

Easy installation

The Direct-Q® range of systems requires no special installation. You can easily set the system up yourself: just connect the system to your tap water supply, plug it in, and insert the SmartPak® purification cartridges. Then, if your configuration includes the [Remote dispenser](#) or a 30- or 60-liter reservoir, follow the simple setup procedures – and your system is ready to use!

Optimized lab space

The compact, all-in-one Direct-Q® water purification unit lets you locate it nearly anywhere in your laboratory, either on or under the benchtop, or wall-installed.*

*Direct-Q® 5 and Direct-Q® 8 systems are designed for use with a 30 l or 60 l external reservoir.



Flow rates to match your requirements

Choose the solution that best meets your lab's requirements with a Direct-Q® system that provides 3, 5, or 8 liters of pure water per hour, and more than 0.5 l of ultrapure water per minute (> 30 l per hour). When not in use, the Direct-Q® system automatically recirculates water in order to maintain high water quality, so that you do not have to wait when sourcing ultrapure water. For your convenience, and to save time, you can also set the system to automatically deliver your selected volume of ultrapure water on demand.

Choice of storage volumes

With its built-in reservoir, the Direct-Q® 3 system can store 6 l of reverse osmosis (RO) water, while Direct-Q® 5 and 8 models are designed for use with a 30- or 60-liter reservoir. Just choose the storage volume that best meets your daily water volume needs.



Fine-tune your water quality

Both pure and ultrapure water

The Direct-Q® range of systems provides you with a convenient and flexible solution for your pure and ultrapure water needs, directly from potable tap water. You'll have access to ultrapure water for your critical applications, and pure water for less critical applications such as general glassware washing or final rinsing – from the same system!

The high quality ultrapure water produced by Direct-Q® systems is suitable for applications such as production of mobile phase for chromatographic separations; preparation of blanks and standard solutions for spectrophotometry, spectroscopy or other analytical techniques; and preparation of buffers for biochemical experiments.

Organic-sensitive applications

Direct-Q® systems are also available with a built-in 185/254 nm UV lamp to reduce the level of organics for critical applications. Water with low TOC provides important benefits to HPLC users such as higher sensitivity and longer column lifetime. The same UV lamp also destroys bacteria.

Application Pak point-of-use polishers

EMD Millipore's range of Application Pak polishers makes it possible to fine-tune your ultrapure water quality to match your research needs. Are your applications sensitive to bacteria, particulates, pyrogens, nucleases, endocrine disruptors or Volatile Organic Compounds? If so, just choose the appropriate final polisher from our range of Application Paks to provide optimal water quality for your requirements.

Please see www.millipore.com/labwater for more information.



Stay focused on your work

Versatile remote dispenser

Designed to fit perfectly into your lab environment, the versatile **Remote dispenser** can be placed up to two meters from your Direct-Q® water purification unit. Select the free-standing or wall-installed model according to your needs – their ergonomony will make either one a welcome addition to your lab, giving you the freedom to focus on your research, while dispensing ultrapure water exactly where you need it. Alternatively, Direct-Q® systems are also available with an integrated dispenser for use on the benchtop.

The design of the bench-installed Direct-Q® systems is also adapted to the height and shape of common laboratory glassware.



EMD Millipore offers more than water

Just the information you need

The intuitive color graphic display shows key system parameters at a glance, enabling easy water quality and maintenance warnings monitoring; the screen rotates for easy viewing wherever the system is located. A short *Quick Reference Guide* is conveniently located in the base of the system for immediate information; the complete *User Manual* can be stored at the back of the water purification unit.

User-friendly maintenance

The all-in-one SmartPak® purification cartridges are easily replaced in just a few minutes.

Watercare Pact Service portfolio

To optimize the performance and lifetime of your water purification system, EMD Millipore offers a complete portfolio of Service plans ranging from a single annual checkup to a full system cover. For more information, please check with your EMD Millipore applications specialist or visit our website: www.millipore.com/labwater



Specifications

Ultrapure (Type I) Product Water Quality*

Resistivity

Production flow rate Direct-Q® 3

Production flow rate Direct-Q® 5

Production flow rate Direct-Q® 8

Instant flow rate
(with Application Pak final filter)

TOC (w/o 185/254 nm UV lamp)

TOC (with 185/254 nm UV lamp)

Particulates (size > 0.22 µm)**

Bacteria**

Endotoxin***(pyrogens)

RNases***

DNases***

* In regular operating conditions

** With Millipak® Express 20 (0.22 µm) membrane filter or with BioPak® ultrafiltration cartridge as final polisher

*** Only with BioPak® ultrafiltration cartridge as final polisher

Direct-Q® Systems

18.2 MΩ·cm @ 25 °C

3 l/h @ 25 °C +/- 15 %

5 l/h @ 25 °C +/- 15 %

8 l/h @ 25 °C +/- 15 %

> 0.5 l/min

< 10 ppb

< 5 ppb

< 1 particulate/ml

< 0.1 cfu/ml

< 0.001 EU/ml

< 0.01 ng/ml

< 4 pg/µl

Pure (Type III) Product Water Quality*

Ionic rejection

Organic rejection for MW > 200

Bacteria and particulates

* In regular operating conditions

> 96 %


> 99 %

> 99 %



System Information	
Dimensions (H x W x D)	54 x 29 x 38 cm (21.3 x 11.4 x 15 in)
Net weight (Direct-Q® 3 system w/o 185/254 nm UV lamp)	8.1 kg (17.9 lb)
Net weight (Direct-Q® 3 system with 185/254 nm UV lamp)	8.6 kg (19.0 lb)
Net weight (Direct-Q® 5, 8 systems with 185/254 nm UV lamp)	7.6 kg (16.7 lb)
Operating weight (Direct-Q® 3 system w/o 185/254 nm UV lamp)	17.6 kg (38.8 lb)
Operating weight (Direct-Q® 3 system with 185/254 nm UV lamp)	18.2 kg (40.1 lb)
Operating weight (Direct-Q® 5, 8 systems with 185/254 nm UV lamp)	12.2 kg (26.8 lb)
Net weight (Remote dispenser)	2.15 kg (4.8 lb)
Operating weight (Remote dispenser)	2.68 kg (5.91 lb)
Built-in reservoir volume	6 l
Electrical feed voltage	100-250 V +/- 10 %
Electrical feed frequency	50-60 Hz +/- 10 %
Tap (feed) water connection	½" Gaz M
Tap (feed) water pressure	0.5 to 6 bar

Available System Configurations

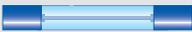Available Configurations	Water Purification Systems		
	Direct-Q® 3	Direct-Q® 5	Direct-Q® 8
UV 185/254 nm 	With/Without	With	With
Remote dispenser 	With/Without	With/Without	With/Without
Built-in 6 l reservoir 	With	Without	Without
30/60 l reservoir* 	Option	Required	Required

* 30 l or 60 l depending on daily water usage



For more information, please visit our website:

www.millipore.com/directq358

Millipore, Direct-Q®, SmartPak®, Millipak® and BioPak® are registered trademarks of Merck KGaA, Darmstadt, Germany.
EMD Millipore and the M mark are trademarks of Merck KGaA.

Lit. No. PB2040ENUS

© 2012 EMD Millipore Corporation, Billerica, MA, U.S.A. All rights reserved.